

DRESSING MATERIAL



MISCELLANEOUS



The Ultimate In Surgical Disposables



Plot No. 10/B, Shyam Ujjawal Ind.Estate, Opp. SBI Bank,
Phase-1, GIDC Vatva, Ahmedabad-382445, Gujarat, India.
Contact No. +91-79-29701333 | Customer Care No. +91-958196169
Email: apexmedical@live.in | dranil@apexmedicaldevices.com
Website: www.apexmedicaldevices.com

ISO 9001:2016 Company | ISO 13485:2016 Certified | CE Certified Company

INFUSION & TRANSFUSION




UROLOGY



ANESTHESIA



GASTROENTEROLOGY

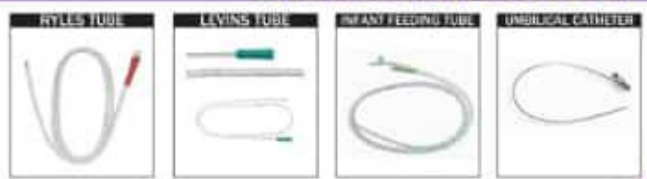

APEX300™

Be Sure - Be Safe

INFUSION SET

- ▶ AS PER IS 12655 ▶ GRAVITY FEED ONLY
- ▶ DISPOSABLE ▶ NON-TOXIC ▶ PYROGEN FREE
- ▶ 20 DROPS = 1ML ± 0.1ML (D.W.)

CAUTION:

- ▶ DISCARD AFTER SINGLE USE.
- ▶ DO NOT USE IF PACKAGE IS OPENED OR DAMAGED.
- ▶ DO NOT STORE AT EXTREME TEMPERATURE AND HUMIDITY.
- ▶ DO NOT USE IF PROTECTIVE CAPS ARE LOOSE OR MISSING.
- ▶ DO NOT USE FOR BLOOD & BLOOD COMPONENTS.

Mfg. Lic No.: G/28/1400

Manufactured By:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC
Vatva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333,
Email: apexmedical@live.in,
web: www.apexmedicaldevices.com

STERILE EO

Quantity: 1 Unit



INSTRUCTION FOR USE:

- ▶ CLOSE THE ROLLER CLAMP
- ▶ REMOVE PROTECTIVE CAP FROM PLASTIC SPIKE & INSERT THE SPIKE AT THE TOP OF THE BOTTLE TO ITS FULL LENGTH
- ▶ SQUEEZE THE DRIP CHAMBER UNTIL IT IS HALF FILLED.
- ▶ REMOVE NEEDLE PROTECTOR HOLDING NEEDLE UPRIGHT. OPEN ROLLER CLAMP AND ALLOW SOLUTION TO PASS UNTIL ALL AIR BUBBLES IN THE TUBE ARE REMOVED. CLOSE THE ROLLER CLAMP.
- ▶ PERFORM VEIN-PUNCTURE AND REGULATE FLOW BY OPENING CLAMP GRADUALLY TO ACHIEVE THE DESIRED FLOW RATE.

LOT	BATCH NO.	: A-2009
MFG/STR. DATE	: 03/2020	
EXP. DATE	: 02/2023	
MAX. R.P. ₹.	: 110=00	(Incl. of all Taxes)

APEX300™

Be Sure - Be Safe

INFUSION SET

- ▶ AS PER IS 12655 ▶ GRAVITY FEED ONLY
- ▶ DISPOSABLE ▶ NON-TOXIC ▶ PYROGEN FREE
- ▶ 20 DROPS = 1ML ± 0.1ML (D.W.)

CAUTION:

- ▶ DISCARD AFTER SINGLE USE.
- ▶ DO NOT USE IF PACKAGE IS OPENED OR DAMAGED.
- ▶ DO NOT STORE AT EXTREME TEMPERATURE AND HUMIDITY.
- ▶ DO NOT USE IF PROTECTIVE CAPS ARE LOOSE OR MISSING.
- ▶ DO NOT USE FOR BLOOD & BLOOD COMPONENTS.

Mfg. Lic No.: G/28/1400

Manufactured By:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC
Vatva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333,
Email: apexmedical@live.in,
web: www.apexmedicaldevices.com

STERILE EO

Quantity: 1 Unit




INSTRUCTION FOR USE:

- ▶ CLOSE THE ROLLER CLAMP
- ▶ REMOVE PROTECTIVE CAP FROM PLASTIC SPIKE & INSERT THE SPIKE AT THE TOP OF THE BOTTLE TO ITS FULL LENGTH
- ▶ SQUEEZE THE DRIP CHAMBER UNTIL IT IS HALF FILLED.
- ▶ REMOVE NEEDLE PROTECTOR HOLDING NEEDLE UPRIGHT. OPEN ROLLER CLAMP AND ALLOW SOLUTION TO PASS UNTIL ALL AIR BUBBLES IN THE TUBE ARE REMOVED. CLOSE THE ROLLER CLAMP.
- ▶ PERFORM VEIN-PUNCTURE AND REGULATE FLOW BY OPENING CLAMP GRADUALLY TO ACHIEVE THE DESIRED FLOW RATE.

LOT	BATCH NO.	: A-2009
MFG/STR. DATE	: 03/2020	
EXP. DATE	: 02/2023	
MAX. R.P. ₹.	: 110=00	(Incl. of all Taxes)

APEX300™ INFUSION SET N.V.

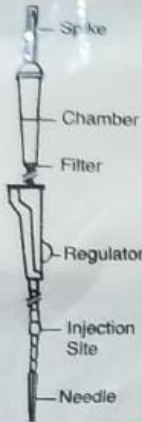

STERILE E.O.

CAUTION

- Do not use if package is open or damaged.
- Do not store in direct sun light at extreme temperature or in high humidity.
- Discard after single use.
- Do not use with blood or blood components.
- Not to be re use.

• AS PER IS 12655 • 20 DROPS = 1ML ± 0.1ML (D.W.)

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC
Vatva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333,
Email : apexmedical@live.in,
web:www.apexmedicaldevices.com



ISO 13485 : 2016

ISO 9001: 2008

- NON PYROGENIC
- GRAVITY FEED ONLY
- NON TOXIC

1 SET

INSTRUCTION FOR USE :

- ▶ CLOSE THE ROLLER CLAMP.
- ▶ REMOVE PROTECTIVE CAP FROM PLASTIC SPIKE & INSERT THE SPIKE AT THE TOP OF THE BOTTLE TO ITS FULL LENGTH.
- ▶ SQUEEZE THE DRIP CHAMBER UNTIL IT IS HALF FILLED.
- ▶ REMOVE NEEDLE PROTECTOR HOLDING NEEDLE UPRIGHT. OPEN ROLLER CLAMP AND ALLOW SOLUTION TO PASS UNTIL ALL AIR BUBBLES IN THE TUBE ARE REMOVED CLOSE THE ROLLER CLAMP.
- ▶ PERFORM VEIN-PUNCTURE AND REGULATE FLOW BY OPENING CLAMP GRADUALLY TO ACHIEVE THE DESIRED FLOW RATE.



Mfg. Lic No.: G/28/1400

APEX300™ INFUSION SET N.V.



STERILE E.O.

CAUTION

- Do not use if package is open or damaged.
- Do not store in direct sun light at extreme temperature or in high humidity.
- Discard after single use.
- Do not use with blood or blood components.
- Not to be re use.

• AS PER IS 12655 • 20 DROPS = 1ML ± 0.1ML (D.W.)

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC
Vatva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333,
Email : apexmedical@live.in,
web:www.apexmedicaldevices.com



ISO 13485 : 2016

ISO 9001: 2008

- NON PYROGENIC
- GRAVITY FEED ONLY
- NON TOXIC

1 SET

INSTRUCTION FOR USE :

- ▶ CLOSE THE ROLLER CLAMP.
- ▶ REMOVE PROTECTIVE CAP FROM PLASTIC SPIKE & INSERT THE SPIKE AT THE TOP OF THE BOTTLE TO ITS FULL LENGTH.
- ▶ SQUEEZE THE DRIP CHAMBER UNTIL IT IS HALF FILLED.
- ▶ REMOVE NEEDLE PROTECTOR HOLDING NEEDLE UPRIGHT. OPEN ROLLER CLAMP AND ALLOW SOLUTION TO PASS UNTIL ALL AIR BUBBLES IN THE TUBE ARE REMOVED CLOSE THE ROLLER CLAMP.
- ▶ PERFORM VEIN-PUNCTURE AND REGULATE FLOW BY OPENING CLAMP GRADUALLY TO ACHIEVE THE DESIRED FLOW RATE.



Mfg. Lic No.: G/28/1400

LOT: BATCH NO. :
MFG./STR. DATE :
EXP. DATE :
MAX. R.P. ₹. :
(Incl. of all Taxes)

IV-68
05/2019
04/2022
185-00

ISO 13485:2016

APPEAL MEDICAL DEVICES
100, 1st Floor, 1st Phase,
Chennai - 600 008, INDIA
Vatva,
Customer:
Email: sales@appealdevices.com
web: www.appealdevices.com

Needle
Site
Injection

LOT: BATCH NO. :
MFG./STR. DATE :
EXP. DATE :
MAX. R.P. ₹. :
(Incl. of all Taxes)

STERILE EO

120 Pouch
Dr. Gauri - 981-10-1310333
medical@live.in
www.appealdevices.com



CE CERTIFIED COMPANY



100 GLOVES
BY WEIGHT

MADE IN MALAYSIA

APEX 300™ INFUSION SET N.V.

- NON PYROGENIC
- NON TOXIC
- GRAVITY FEED ONLY
- Needle Size 21 G



STERILE E.O.

CAUTION

- Do not use if package is open or damaged.
- Do not store in direct sun light at extreme temperature or in high humidity.
- Discard after single use.
- Do not use with blood or blood components.
- Not to be re used.

• AS PER IS 12855 • 20 DROPS = 1ML ± 0.1ML (D.W.)

INSTRUCTION FOR USE

- Close the roller clamp.
- Remove protector from spike & insert the spike at the top of the bottle to its full length.
- Squeeze the drip chamber until it is fluid half.
- Remove needle protector holding needle upright.
- Open roller clamp & allow solution to pass until all air bubbles in the tube are removed.
- Close the roller clamp.
- Perform Vein puncture and regulate flow by opening flow regulator.
- Mfg. Lic No.: G/28/1400

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujwala Ind Park,
Opp. State Bank of India, Phase-1, GIDC
Valva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333.
Email: apexmedical@live.in.
web: www.apexmedicaldevices.com

CE

ISO 13485: 2016

ISO 9001: 2008



1 SET

BATCH NO. : IV-014
INFO. STR. DATE : 02/2020
EXP. DATE : 01/2023
MAX. R.P. : 149=00
(Incl. of all Taxes)

APEX 300™ INFUSION SET N.V.

- NON PYROGENIC
- NON TOXIC
- GRAVITY FEED ONLY
- Needle Size 21 G



STERILE E.O.

CAUTION

- Do not use if package is open or damaged.
- Do not store in direct sun light at extreme temperature or in high humidity.
- Discard after single use.
- Do not use with blood or blood components.
- Not to be re used.

• AS PER IS 12855 • 20 DROPS = 1ML ± 0.1ML (D.W.)

INSTRUCTION FOR USE

- Close the roller clamp.
- Remove protector from spike & insert the spike at the top of the bottle to its full length.
- Squeeze the drip chamber until it is fluid half.
- Remove needle protector holding needle upright.
- Open roller clamp & allow solution to pass until all air bubbles in the tube are removed.
- Close the roller clamp.
- Perform Vein puncture and regulate flow by opening flow regulator.
- Mfg. Lic No.: G/28/1400

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujwala Ind Park,
Opp. State Bank of India, Phase-1, GIDC
Valva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333.
Email: apexmedical@live.in.
web: www.apexmedicaldevices.com

CE

ISO 13485: 2016

ISO 9001: 2008



1 SET

BATCH NO. : IV-014
INFO. STR. DATE : 02/2020
EXP. DATE : 01/2023
MAX. R.P. : 149=00
(Incl. of all Taxes)





APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES



THE ULTIMATE IN SURGICAL DISPOSABLE PRODUCTS

APEX
MEDICAL DEVICES



CE CERTIFIED COMPANY

APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES


APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES

THE ULTIMATE IN SURGICAL DISPOSABLE PRODUCTS

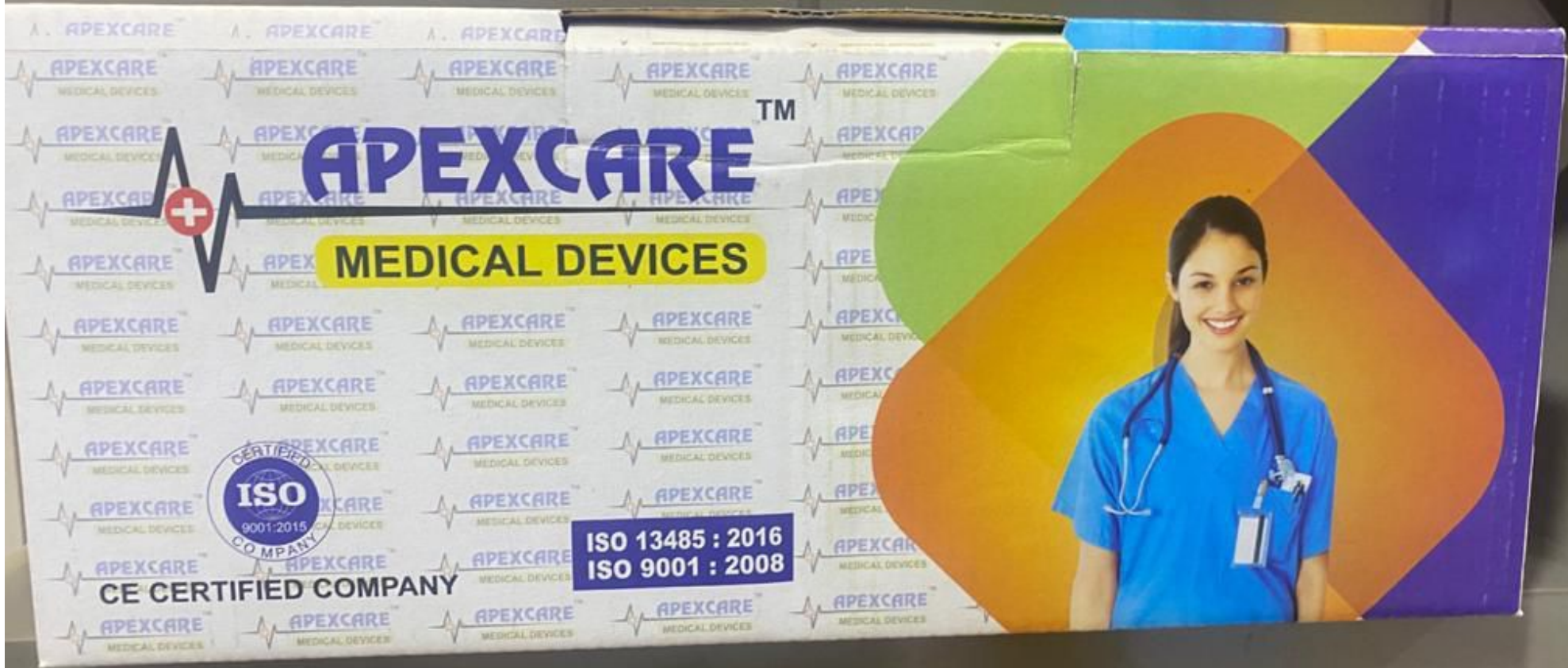

APEXCARETM

MEDICAL DEVICES

CERTIFIED
ISO
9001:2015
COMPANY

CE CERTIFIED COMPANY

ISO 13485 : 2016
ISO 9001 : 2008

INFUSION SET

STERILE	☒
NON-TOXIC	☒
PYROGEN FREE	☒
SINGLE USE ONLY	☒
GRAVITY FEED ONLY	☒

Qty.: 50 Units



CAUTION

- Do not use if primary packaging is damaged.
- Discard after single use.
- Do not re-sterilize.
- Store in cool & dry place.

Mfg.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind. Park,
Opp. State Bank of India, Phase-1, GIDC
Valva, Ahmedabad-382445, Gujarat, India.
Customer Care : +91-79-29701333,
E-mail : apexmedical@live.in,
Website : www.apexmedicaldevices.com

00300



APEX MEDICAL DEVICES

VACOVEIN INFUSION SET

- NON PYROGENIC
- NON TOXIC
- GRAVITY FEED ONLY
- Needle Size 21 G



STERILE E.O.

CAUTION

- Do not use if package is open or damaged.
- Do not store in direct sun light at extreme temperature or in high humidity.
- Discard after single use.
- Do not use with blood or blood components.
- Not to be re use.
- AS PER IS 12655
- 20 DROPS = 1ML $\pm$ 0.1ML (D.W.)

INSTRUCTION FOR USE

- Close the roller clamp.
- Remove protector from spike & insert the spike at the top of the bottle to its full length.
- Squeeze the drip chamber until it is fluid half.
- Remove needle protector holding needle upright.
- Open roller clamp & Air vent snap cap to remove air from tubing.
- Close the roller clamp.
- Perform Vein puncture and regulate flow by opening flow regulator.

Mfg. Lic No.: G/28/1400

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC-Valva,
Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333.

WITH AIRVENT



CE
ISO 13485: 2016
ISO 9001: 2008

1 SET

BATCH NO. 1U-001
MFG/STR. DATE 01/2021
EXP. DATE 12/2023
MAX. R.P. ₹. 155/-
(Incl. of all Taxes)

VACOVEIN INFUSION SET

- NON PYROGENIC
- NON TOXIC
- GRAVITY FEED ONLY
- Needle Size 21 G



STERILE E.O.

CAUTION

- Do not use if package is open or damaged.
- Do not store in direct sun light at extreme temperature or in high humidity.
- Discard after single use.
- Do not use with blood or blood components.
- Not to be re use.
- AS PER IS 12655
- 20 DROPS = 1ML $\pm$ 0.1ML (D.W.)

INSTRUCTION FOR USE

- Close the roller clamp.
- Remove protector from spike & insert the spike at the top of the bottle to its full length.
- Squeeze the drip chamber until it is fluid half.
- Remove needle protector holding needle upright.
- Open roller clamp & Air vent snap cap to remove air from tubing.
- Close the roller clamp.
- Perform Vein puncture and regulate flow by opening flow regulator.


Mfg. Lic No.: G/28/1400

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC-Valva,
Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333.


WITH AIRVENT

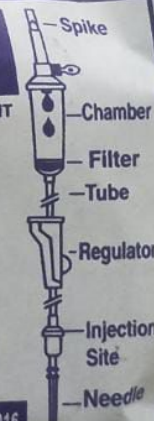

CE
ISO 13485: 2016
ISO 9001: 2008

1 SET

BATCH NO. 1U-001
MFG/STR. DATE 01/2021
EXP. DATE 12/2023
MAX. R.P. ₹. 155/-
(Incl. of all Taxes)

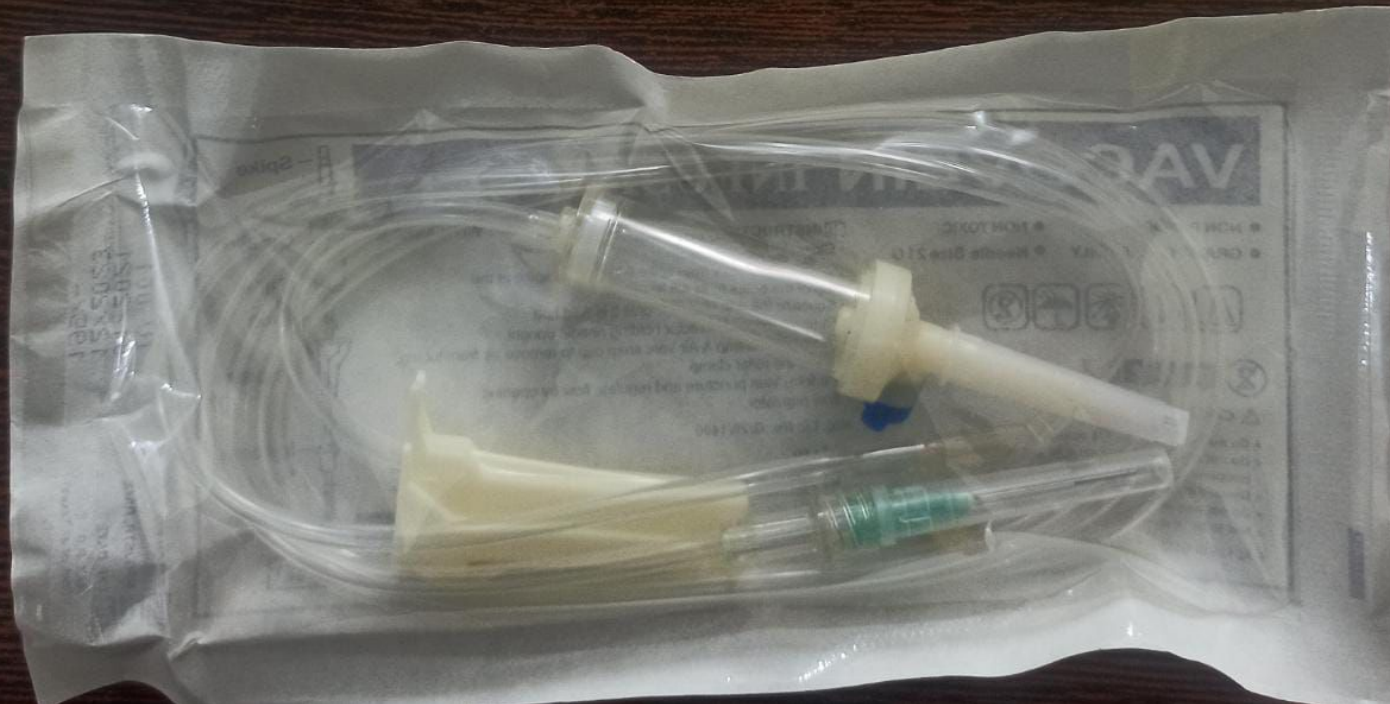

VACOFIX

MICRODIP SET

- NON PYROGENIC
- NON TOXIC
- GRAVITY FEED ONLY
- Needle Size 23 G



 **STERILE E.O.**

CAUTION :

- This is for single use discard after use.
- Do not store in direct sunlight at extreme temperature or in high humidity.
- Do not use if package is open or damaged.
- Do not use with blood or blood components.
- Mfg. Lic No.:G/28/1400

• AS PER IS 12655 • 60 DROPS = 1ML $\pm$ 0.1ML (D.W.)

INSTRUCTION FOR USE :

1. Remove the outer cap of the container and disinfect top.
2. Firmly close regulator.
3. Remove spike cover and insert the spike at the top of the container to its full length, turn the container upside down and hang it at suitable height.
4. Squeeze the dripchamber until it is filled upto the mark.
This prevents air bubbles from entering the fluid path.
5. Remove the protector of the needle Gently release the regulator in order to let the solution flow, Keeping the needle above the container level. Gradually lower the needle in such a way that no air bubbles are left in the tube. Let a few drops of solution pass through.
6. Close the regulator.
7. Perform veinpuncture and release the regulator to achieve the desired flow rate.

Mfr.

APAR MEDICAL DEVICES

PLOT No. 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC-Vatva,
Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333.

CE

ISO 18485 : 2016

ISO 9001: 2008

1 SET



MATCH DATE
MFG/STR. DATE : 08/2020
EXP. DATE : 07/2023
MAX. R.P. : 195-00
(Incl. of all Taxes)

- Store in a cool & dry place
- Discard the set after single use
- Do not use if protective caps are loose or missing
- Do not use if package is opened or damaged

CAUTION:

STERILE EO



AS PER IS 8854

Needle size : 18g

50 Drops per ml ± 0.1 ml

NON TOXIC, PYROGEN FREE, SINGLE USE

APEX 300 BLOOD ADMINISTRATION SET

Customer Care: +91-38-38301333
 Ahmedabad-383442, Gujarat, India
 Opp. State Bank of India, Phase-1, CIDC-Vastra,
 Plot No 101B, Shyam Ullas Ind Park

APEX MEDICAL DEVICES

Model No.: Apex 300

1. Clamp gradually to achieve the desired flow rate
2. Perform vein-puncture and regulate flow by opening bubbles in the tube are removed close the roller clamp
3. Open roller clamp & allow blood to pass until all air
4. Remove needle protector holding needle upright
5. Squeeze the drip chamber until it is half filled
6. Remove Protective cap from plastic spike & insert the
7. Close the roller clamp

INSTRUCTIONS FOR USE



120 3004: 3008
 120 13482: 3016



APEX 300TM BLOOD ADMINISTRATION SET

NON TOXIC, PYROGEN FREE, SINGLE USE
20 Drops per ml. ± 0.1 ml.

Needle Size : 18G
AS PER IS 9824



STERILE E.O.

CAUTION:

- Do not use if package is opened or damaged
- Do not use if protective caps are loose or missing
- Discard the set after single use
- Store in a cool & dry place

INSTRUCTIONS FOR USE


1. Close the roller clamp
2. Remove Protective cap from plastic spike & insert the spike at the top of the bottle to its full length
3. Squeeze the drip chamber until it is half filled
4. Remove needle protector holding needle upright
5. Open roller clamp & allow blood to pass until all air bubbles in the tube are removed close the roller clamp
6. Perform vein-puncture and regulate flow by opening clamp gradually to achieve the desired flow rate

Mfg. Lic No.: G/28/1400

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC-Valva,
Ahmedabad-382445, Gujarat, India
Customer Care: +91-79-29701333

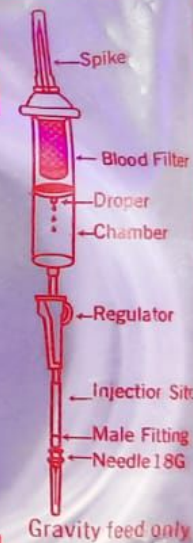

1 SET

LOT: 42152
MFG/STR. DATE: 06/2024
EXP. DATE: 05/2024
MAX. R.P. 160-00
(Incl. of all Taxes)



APEX 300™

MEASURED VOLUME SET

(For Pediatric Use Only)

STERILE, NON TOXIC, PYROGEN FREE, SINGLE USE ONLY
"Gravity Feed Only"

SPECIAL FEATURES :

- * 60 Drops equivalent to 1 ml $\pm$ 0.1 ml (1g $\pm$ 0.1g)
- * Floating Shut-Off Valve.

CAUTION :

- DO NOT USE IF PACKING IS OPEN OR DAMAGED.
- DO NOT USE IF PROTECTIVE CAPS ARE LOOSE OR MISSING.
- DO NOT USE FOR BLOOD AND BLOOD COMPONENTS.
- USE RUBBER PARTS ONLY FOR ADDITIONAL DRUG ADMINISTRATION.
- DESTROY THE INFUSION SET AFTER USE.

and Mr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC
Valva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333,
Email: apexmedical@live.in,
web: www.apexmedicaldevices.com

CE

ISO 13485 : 2016
ISO 9001: 2008

INSTRUCTIONS FOR USE :


1. Close all controllers.
2. Remove spike protector and insert firmly into the closing stopper of the infusion container.
3. Suspend infusion container with attached set. Open Upper clamp and allow approximately 30 ml. solution to flow into the graduated chamber. Close the upper clamp.
4. Gently squeeze and release drip chamber until it is approximately one third full.
5. Remove vein needle protector. Slightly open flow regulator to clear air from tubing and the vein needle. Close the flow regulator.
6. Open the upper clamp & allow solution to flow into graduated chamber until the desired volume is obtained. Close the upper clamp. Perform vein puncture.
7. Gradually open the flow regulator. Adjust drop rate & control the solution flow with the flow regulator.
8. Hang the complete set with the help of the hanger on the other side of the I.V. stand, to prevent accidental fall of the Infusion Set.
9. The floating valve shall shut-off the flow when solution level in the graduated chamber comes to zero.
10. When more solution is required close flow regulator & open upper clamp to fill graduated chamber upto the desired level. Close the upper clamp and squeeze the drip chamber gently to float the shut-off valve and then restart infusion as in step - 8.
11. Close air vent Caps during periods of interrupted infusion therapy.

STERILE EO



Mfg. Lic No.: G/28/1400

LOT: 25
A-2172
MFG. TR. DATE: 09/2021
EXP. DATE: 08/2022
MAX. P.P. 210/-
(per set of 100 sets)

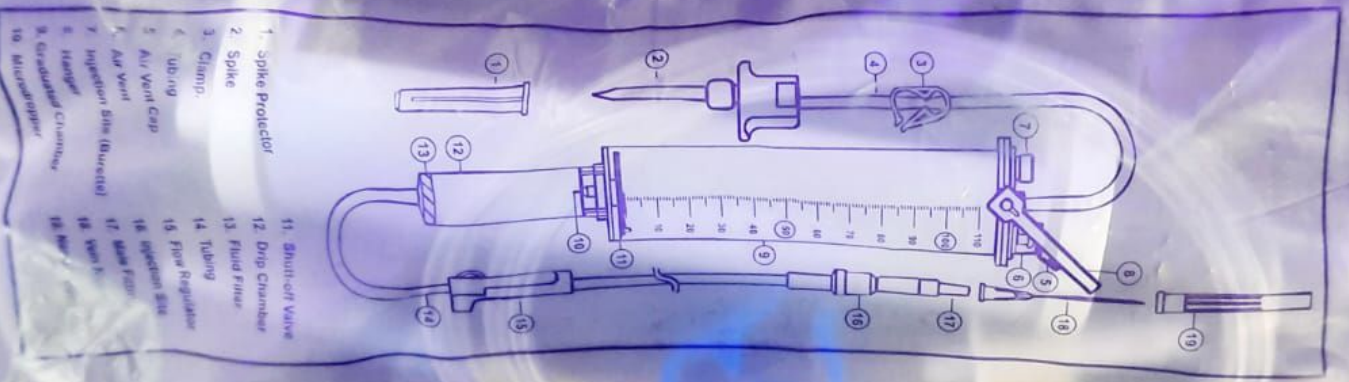

- 1. Spike Protector
- 2. Spike
- 3. Clamp
- 4. Tubing
- 5. Air Vent Cap
- 6. Air Vent
- 7. Injection Site (Burette)
- 8. Hanger
- 9. Graduated Chamber
- 10. Microdrip
- 11. Shut-off Valve
- 12. Drip Chamber
- 13. Fluid Filter
- 14. Tubing
- 15. Flow Regulator
- 16. Injection Site
- 17. Main Flow
- 18. Vent
- 19. Microdrip

INFANT MUCUS EXTRACTOR

CAUTIONS :

- THE SET USE ONCE ONLY
- ENSURE PACKAGE IS UNDAMAGED IF FOUND DISCARD THE SET.
- USE UNDER EXPERT MEDICAL SUPERVISION ONLY
- KEEP IN COOL DRY & DARK PLACE

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC-Vatva,
Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333.
Mfg. by : P.M. Medisurge, Dinubhai Complex,
GIDC Vatva, Ahmedabad-382445, Gujarat, INDIA



ISO 13485 : 2016

ISO 9001: 2008



STERILE EO

- Non-Toxic
- Non Pyrogenic
- E.O. Gas Sterilised
- Disposable

Mfg. Lie. No.: G/28/1374

1 UNIT



LOT	BATCH NO.	:	10
MFG/STR. DATE	:	12/2019	
EXP. DATE	:	11/2022	
MAX. R.P. ₹.	:	75/-	
(Incl. of all Taxes)			

MAX. R.P.RS. 22=00
(INCL. OF ALL TAXES)

LOT BATCH NO. : EX-002
MFG/STR. DATE : 08/2020
EXP. DATE : 07/2023

Tear Here

APEX CARE

UMBILICAL CORD CLAMP

CAUTION

- Disposable • Single Use
- Do not store at extreme temperature and at high humidity.
- Do not reuse. Reuse of the devices may cause injury, illness or even death because of biological contamination.



1 UNIT



ISO 13485 : 2016
ISO 9001 : 2008

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC-Vatva,
Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333,


Customer Care: +91-79-29701333,
Ahmedabad-382445, Gujarat, India.
Opp. State Bank of India, Phase-1, GIDC-Vatva,
Plot No 10/B, Shyam Ujjawal Ind Park,

APEX MEDICAL DEVICES

Mfr.:



ISO 13485 : 2016
ISO 9001 : 2008

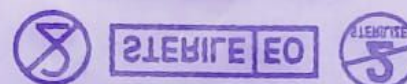

1 UNIT

- illness or even death because of biological contamination.
- Do not reuse. Reuse of the devices may cause injury, illness or even death because of biological contamination.
 - Do not store at extreme temperature and at high humidity.
 - Disposable • Single Use

CAUTION



EXP. DATE : 07/2023
MFG/STR. DATE : 08/2020

MAX. R.P.RS. 22=00
(INCL. OF ALL TAXES)



ISO 13485 : 2016
ISO 9001 : 2015
CE CERTIFIED COMPANY

Foley Balloon Catheter



Foley Balloon Catheter



ISO 13485 : 2016
ISO 9001 : 2015
CE CERTIFIED COMPANY

Foley Balloon Catheter



DIRECTIONS:
Following aseptic catheterization inflate the retention balloon by inserting the syringe to tightly into the valve and introduce the appropriate volume of sterile water.
Recommended Inflation Capacities

WARNING:
Do not use ointments or lubricants having a petroleum base.
Do not re-use. Reuse of the device may cause injury, illness or even death because of biological contamination.

Balloon Size	Minimum Volume	Maximum Volume
3-5 ml	2 ml	5 ml
5-15 ml	5 ml	15 ml
15-30 ml	15 ml	30 ml
30-50 ml	30 ml	50 ml



CAUTION

- This device to sale by or on the order of a physician
- Avoid freezing and excessive heat.
- Store in a cool and dry place.
- Do not exceed recommended capacity
- Sterility guaranteed if package is undamaged

NATURAL RUBBER LATEX SILICON ELASTOMER COATED

MANUFACTURED BY:
APEX MEDICAL DEVICES
Plot - 105, Shree Jagad Industrial Park,
Opp. S.S.I. Zone, Phase - 1, S.I.D.C. Sahas,
Ahmednagar 432445, Gujarat, India.
E-mail: apexmedical@rediffmail.com
apexmedical@apexmedicaldevices.com
Web: www.apexmedicaldevices.com

TYPE:
SIZE (F/GFR):
BATCH NO.:
MFG DATE:
EXP DATE:
Max M.R.P. ₹:



ISO 13485 : 2016
ISO 9001 : 2015
CE CERTIFIED COMPANY

Foley Balloon Catheter



Foley Balloon Catheter



ISO 13485 : 2016
ISO 9001 : 2015
CE CERTIFIED COMPANY

Foley Balloon Catheter



CAUTION

- This device to sale by or on the order of a physician
- Avoid freezing and excessive heat.
- Store in a cool and dry place.
- Do not exceed recommended capacity
- Sterility guaranteed if package is undamaged

NATURAL RUBBER LATEX SILICON ELASTOMER COATED



APEXCARE™
Be Sure - Be Safe

URINE COLLECTION BAG

ECONOMY

FEATURES:

- Urine Bag Made from clinical grade soft PVC, semi transparent bags.
- 2 Ltr. Capacity
- A Non return valve has been provided in the Urine Bag which prevents back flow, thus offering added protection.
- E.T.O. Sterilized.
- Biologically tested and ready for single use.
- Extra soft & long tube.
- Easy to fit any size of Foley's catheters,



CAUTION:

- Do not store in direct sunlight at extreme temperature or in high humidity.
- Do not use if package is opened or damaged or protective cap is detached.
- Do not use for more than five days.
- Do not reuse : it may cause infection.



CE ISO 13485 : 2016
ISO 9001 : 2015

LOT BATCH NO. : AUB-005
MFG./STR. DATE : 01/2021
EXP. DATE : 12/2023
MAX. R.P. ₹. : 125.00
(Incl. of all Taxes)

Manufactured By:
APEX MEDICAL DEVICES

Plot No. 12, Sector 10, Gurgaon - 122 002
Dist. Gurgaon, Haryana
Contact No. : 0122-2211111
Fax : 0122-2211112
Email : apexmedical@rediffmail.com
www.apexmedical.com



Urine Drainage Bag
With
Non-Return Valve
DOUBLE SEAL

Wider Bag and Double Scale ensure better reading of diuresis.

To measure small volumes of 25ml to 100 ml hold near outlet corner and allow bag to hang freely.

Full pressure compensation is provided to ensure free entry of urine into bag.

I.B. : Urine must not exceed MAX. CAPACITY level.

DOUBLE SEAL

1500 1000 500 400 300

Patient :
Bed :
Ward :
Doctor :
Date :

APEXCARE™
Be Sure - Be Safe

URINE COLLECTION BAG

ECONOMY

FEATURES:

- Urine Bag Made from clinical grade soft PVC, semi transparent bags.
- 2 Ltr. Capacity
- A Non return valve has been provided in the Urine Bag which prevents back flow, thus offering added protection.
- E.T.O. Sterilized.
- Biologically tested and ready for single use.
- Extra soft & long tube.
- Easy to fit any size of Foley's catheters,



CAUTION:

- Do not store in direct sunlight at extreme temperature or in high humidity.
- Do not use if package is opened or damaged or protective cap is detached.
- Do not use for more than five days.
- Do not reuse, it may cause infection.



CE ISO 13485 : 2016
ISO 9001 : 2015

LOT BATCH NO. : AUB-005
MFG/STR. DATE : 01/2021
EXP. DATE : 12/2023
MAX. R.P. ₹. : 125.00
(Incl. of all Taxes)



APEXCARE™
Be Sure - Be Safe

URINE COLLECTION BAG

ECONOMY

FEATURES:

- Urine Bag Made from clinical grade soft PVC, semi transparent bags.
- 2 Ltr. Capacity
- A Non return valve has been provided in the Urine Bag which prevents back flow, thus offering added protection.
- E.T.O. Sterilized.
- Biologically tested and ready for single use.
- Extra soft & long tube.
- Easy to fit any size of Foley's catheters,




CAUTION:

- Do not store in direct sunlight at extreme temperature or in high humidity.
- Do not use if package is opened or damaged or protective cap is detached.
- Do not use for more than five days.
- Do not reuse, it may cause infection.



CE ISO 13485 : 2016
ISO 9001 : 2015

LOT BATCH NO. : AUB-005
MFG/STR. DATE : 01/2021
EXP. DATE : 12/2023
MAX. R.P. ₹. : 125.00
(Incl. of all Taxes)

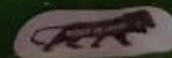

Tear Here

APEXCARE

URINE BAG

FEATURES:

- Urine Bag Made from clinical grade soft PVC, semi transparent bags.
- 2 Ltr. Capacity
- A Non return valve has been provided in the Urine Bag which prevents back flow, thus offering added protection.
- E.T.O. Sterilized.
- Biologically tested and ready for single use.
- Extra soft & long tube.
- Easy to fit any size of Foley's catheters,



CAUTION:

- Do not store in direct sunlight at extreme temperature or in high humidity.
- Do not use if package is opened or damaged or protective cap is detached.
- Do not use for more than five days.
- Do not reuse : it may cause infection.

CE

ISO 13485 : 2016
ISO 9001 : 2015

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC-Vatva,
Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333.

LOT	BATCH NO.	FAUE-002
MFG/STR. DATE	07/2019	
EXP. DATE	06/2022	
MAX. R.P. ₹.	195.00	
(Incl. of all Taxes)		



FEATURES:

- Urine Bag Made from clinical grade soft PVC, semi transparent bags.
- 2 Ltr. Capacity
- A Non return valve has been provided in the Urine Bag which prevents back flow, thus offering added protection.
- E.T.O. Sterilized.
- Biologically tested and ready for single use.
- Extra soft & long tube.
- Easy to fit any size of Foley's catheters.



CAUTION:


- Do not store in direct sunlight or extreme temperature or in high humidity.
- Do not use if package is opened or damaged or protective cap is detached.
- Do not use for more than one time.
- Do not reuse : it may cause infection.

CE

ISO 13485 : 2016
ISO 9001 : 2015

BATCH NO. : AUB-002
MFG. DATE : 07/2019
EXP. DATE : 08/2022
MAX. R.P. : 105.00
(Per. of 10/1000)

APEX MEDICAL DEVICES

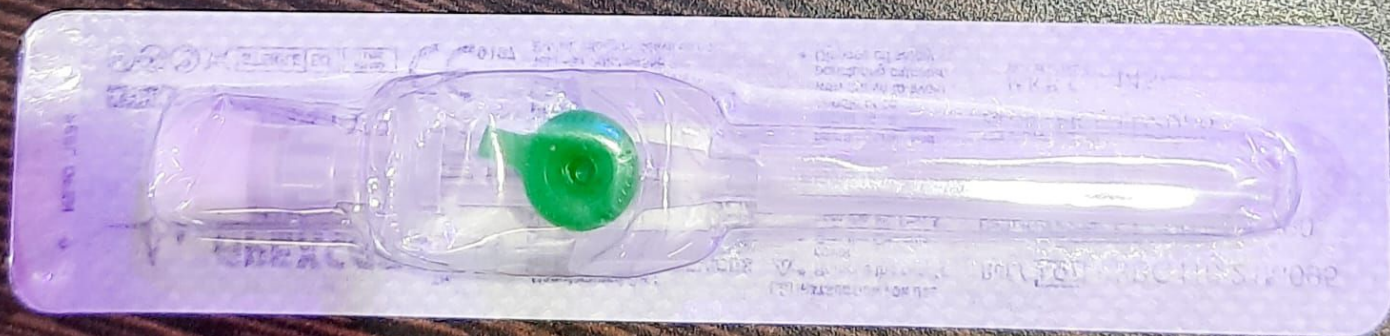

 APEXCARETM
I.V. Cannula
With Wings & Injection Port

RADIO-OPAQUE CATHETER



100 PCS.

STERILE EO

ISO
13485

CE 0197



APEXCARETM

I.V. Cannula

With Wings & Injection Port

WARNING

- The use of this products is restricted a qualified doctor or a paramedic.
- The products should be used according to the instructions for use.
disclaim any responsibility for possible consequences resulting from improper use.
- Store in dry & cool place. Do not expose to heat or direct sunlight.
- The products should not be reprocessed.
- Visually inspect and carefully check the products and packaging before use because Improper transport and handling may cause structural and / or functional damage to device or packaging.
- The products is non -toxic, & non-pyrogenic till the package has not been opened or damaged.
- Do not clean or resterilise the products after opening the packet.
- The product should be used immediately after opening the packaging.
- Perform routine monitoring & venipuncture site maintenance according to medical norms.
- Discard after use.  Mfg. by : SPM Medicare Pvt. Ltd., Noida, Mfg. Lic. No. :15/SC/P of 2017

24 G



013 1 8904379 80263 2 (13) 211200 (17) 281100 (10) 3BG124/21M098 (21) 78398

Manufactured For : **APEX MEDICAL DEVICES**



Plot : 10/B, Shyam Ujjawal Industrial Park,
Opp. S.B.I. Bank, Phase-1, G.I.D.C, Vatva,
Ahmedabad-382445 Gujarat, India.
E-mail : apexmedical@live.in,
domestic@apexmedicaldevices.com,
Web : www.apexmedicaldevices.com
Customer Care No. : +91-9624699959, 9598196169

Ref. / **LOT** : 3BG124/21M098
Cath./G./F (mm) : FEP R/O/24/23
 MFG Dt. : 12/2021
 EXP. Dt. : 11/2026 -
M.R.P. ₹ : **16500/-**
(Incl. of all taxes)





NON TOXIC & NON PYROGENIC
DO NOT USE PACKING IS OPENED OR DAMAGED
SCALP VEIN SET

② **STERILE E.O.**

B.No.SV-47
MFG.07-2019
EXP.06-2022
M.R.P.Rs.26.00
SIZE: 22 G

Mfg for :

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Phase-1,
GIDC Vatva, Ahmedabad-382445, Gujarat,
India. Customer Care: +91-79-29701333,
Mfg by : P.H. Medisurge Vatva, Ahmedabad
Mfg. Lic. No. G/28/1374



NON TOXIC & NON PYROGENIC
DO NOT USE PACKING IS OPENED OR DAMAGED
SCALP VEIN SET

② **STERILE E.O.**

Mfg for :

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Phase-1,
GIDC Vatva, Ahmedabad-382445, Gujarat,
India. Customer Care: +91-79-29701333,
Mfg by : P.H. Medisurge Vatva, Ahmedabad
Mfg. Lic. No. G/28/1374



NON TOXIC & NON PYROGENIC
DO NOT USE PACKING IS OPENED OR DAMAGED
SCALP VEIN SET

② **STERILE E.O.**

B.No.SV-46
MFG.12-2018
EXP.11-2021
M.R.P.Rs.26.00
SIZE: 24 G

Mfg for :

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Phase-1,
GIDC Vatva, Ahmedabad-382445, Gujarat,
India. Customer Care: +91-79-29701333,
Mfg by : P.H. Medisurge Vatva, Ahmedabad
Mfg. Lic. No. G/28/1374



NON TOXIC & NON PYROGENIC
DO NOT USE PACKING IS OPENED OR DAMAGED
SCALP VEIN SET

② **STERILE E.O.**

Mfg for :

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Phase-1,
GIDC Vatva, Ahmedabad-382445, Gujarat,
India. Customer Care: +91-79-29701333,
Mfg by : P.H. Medisurge Vatva, Ahmedabad
Mfg. Lic. No. G/28/1374





ISO 13485 : 2016
CE CERTIFIED COMPANY



100 PCS



SCALP VEIN SET

Short Bevel Siliconized Needle




ISO 13485 : 2016
CE CERTIFIED COMPANY



SCALP VEIN SET

Short Bevel Siliconized Needle

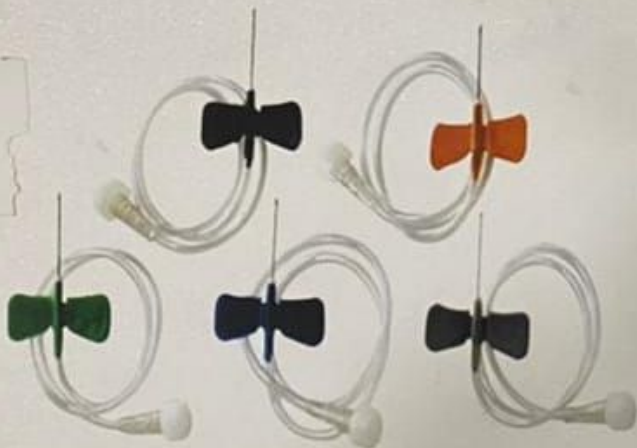

ISO 13485 : 2016
CE CERTIFIED COMPANY



100 PCS



9.1M

MICROPOROUS SURGICAL TAPE

Adhesive Tape Non-Woven Support BP

HYPOALLERGENIC
PAINLESS REMOVER
HIGHLY PERMEABLE TO AIR AND WATER VAPOUR OPTIMUM ADHESION

APEXCARE

Be Sure - Be Safe



ELASTIC ADHESIVE BANDAGE B.P.

APEXFIXTM

I.V.CANNULA FIXATOR

Extra strip for
additional Adhesion



Size : Medium
Batch No. :
Mfg. Date : CF20031
Exp. Date : 3/2020
Max. Retail Price ₹ : 2/2023
(Incl. of all taxes) 70.0
Mfg. Lic. No. : G/25/2104

Mfg.For. :

APEX MEDICAL DEVICES

Plot No.10/B, Shyam Ujjawal Ind. Park, Opp. State Bank of India,
Phase-1, GIDC-Vatva, Ahmedabad-382445, Gujarat, India
Customer Care : +91 79-29701333, Email : apexmedical@live.in
Mfg.By : SteriCott Health Care, Dist: Gandhinagar. (Gujarat)

ELASTIC ADHESIVE BANDAGE B.P.

APEXFIXTM

I.V.CANNULA FIXATOR

Extra strip for
additional Adhesion



Size : Medium
Batch No. : CF20031
Mfg. Date : 3/2020
Exp. Date : 2/2023
Max. Retail Price ₹ : 70.0
(Incl. of all taxes)
Mfg. Lic. No. : G/25/2104

Mfg.For. :

APEX MEDICAL DEVICES

Plot No.10/B, Shyam Ujjawal Ind. Park, Opp. State Bank of India,
Phase-1, GIDC-Vatva, Ahmedabad-382445, Gujarat, India
Customer Care : +91 79-29701333, Email : apexmedical@live.in
Mfg.By : SteriCott Health Care, Dist: Gandhinagar. (Gujarat)

REDMI NOTE 8 PRO
AI QUAD CAMERA

ELASTIC ADHESIVE BANDAGE B.P.

APEXFIXTM

I.V.CANNULA FIXATOR

Extra strip for
additional Adhesion



Size	: Medium
Batch No.	: CF20031
Mfg. Date	: 3/2020
Exp. Date	: 2/2023
Max. Retail Price ₹ (incl. of all taxes)	: 70.0
Mfg. Lic. No.	: G/25/2104

Mfg. For.:

APEX MEDICAL DEVICES

Plot No.10/B, Shyam Ujjawal Ind. Park, Opp. State Bank of India,
Phase-1, GIDC-Vatva, Ahmedabad-382445, Gujarat, India
Customer Care: +91 79-29701333, Email: apexmedical@live.in
Mfg By: SteriCot Health Care, Dist: Gandhinagar (Gujarat)

ELASTIC ADHESIVE BANDAGE B.P.

APEXFIX

I.V.CANNULA FIXATOR

Extra strip for
additional Adhesion




Qty.:
100
Nos.



**STERILE LATEX
SURGICAL GLOVES**



CE CERTIFIED COMPANY

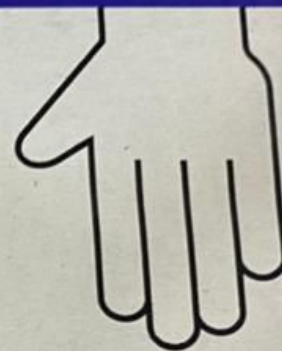

50 PAIRS



OPEN HERE

APEXCARE™

Sterile Latex Surgical Gloves Powdered

FEATURES:

- Pre-powdered with bio absorbable USP grade corn starch.
- Micro textured surface for better instrument grip.
- Anatomically shaped with beaded cuff end.
- Conforms to IS 13422, ASTM D3577 & EN455-1,2,3
- Sterilized by Ethylene Oxide.
- Strictly adhered to GMP guidelines.

CAUTION:

- This product contains natural rubber latex which may cause allergic reactions in some individuals, in case of latex allergy, discontinue use and consult a physician.
- After donning remove powder by wiping gloves thoroughly with a wet sterile sponge, sterile wet towel or other effective methods.
- Sterile until package is open or damaged.
- Store in cool & dry place.
- Do not re-sterilize.



ONE PAIR



FOR SINGLE USE ONLY

STERILE



CE

ISO 13485 : 2016

ISO 9001 : 2015

SIZE: 6.5

LOT

SG005

EXP

05/2021

EXP

04/2024

Max. R.P. ₹ : 125/-
(Incl. of all Taxes)

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Opp. State Bank of India, Phase-1, GIDC
Valva, Ahmedabad-382445, Gujarat, India. Customer Care: +91-79-29701333, 9598196169
Email: apexmedical@live.in, web: www.apexmedicaldevices.com

OPEN HERE

APEXCARE™

Sterile Surgical Latex Gloves Powder-Free

FEATURES:


- Micro textured surface for better instrument grip.
- Anatomically shaped with beaded cuff end.
- Conforms to IS 13422, ASTM D3577 & EN455-1,2,3
- Sterilized by Ethylene Oxide.
- Strictly adhered to GMP guidelines.

CAUTION:

- This product contains natural rubber latex which may cause allergic reactions in some individuals, in case of latex allergy, discontinue use and consult a physician.
- Sterile until package is open or damaged.
- Store in cool & dry place.
- Do not re-sterilize.



ONE PAIR

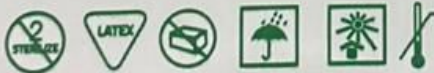

FOR SINGLE USE ONLY
STERILE

CE



ISO 13485 : 2016
ISO 9001: 2015

SIZE : 7.5

LOT SG002
09/2020
08/2023

Max. R.P. ₹ : 125/-
(Incl. of all Taxes)

Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Opp. State Bank of India, Phase-1, GIDC
Valva, Ahmedabad-382445, Gujarat, India. Customer Care: +91-79-29701329
Email : apexmedical@live.in, web: www.apexmedicaldevices.com

APEXCARE™


STERILE SURGICAL RUBBER GLOVES POWDERED

FEATURES:

- Pre-powdered with bio absorbable USP grade corn starch.
- Micro textured Surface for better instrument grip.
- Anatomically shaped with beaded cuff end.
- Conforms to IS 13422, ASTM D3577 & EN455-1,2,3.
- Sterilized by Ethylene Oxide.
- Strictly adhered to GMP guidelines.

CAUTION:

- This product contains natural rubber latex which may cause allergic reactions in some individuals, in case of latex allergy, discontinue use and consult a physician.
- After donning remove powder by wiping gloves thoroughly with a wet sterile sponge, sterile wet towel or effective methods.
- Sterile until package is open or damaged.
- Store in cool & dry place.
- Do not re-sterilize.

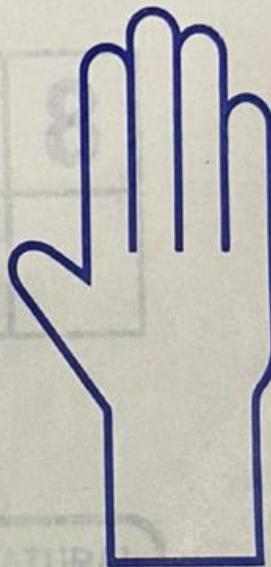

ONE PAIR



FOR SINGLE USE ONLY

STERILE



ISO 13485 : 2016

ISO 9001 : 2015

SIZE : 6.5

LOT : SG 012

EXP : 12-2020

EXP : 11-2023

M.R.P. Rs. : ₹ 125/- Per Pair
(Incl. of all Taxes)



Mfr.:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Opp. State Bank Of India, Phase-1, GIDC
Valva, Ahmedabad-382445, Gujarat, India. Customer Care: +91-79-29701333
Email : apexmedical@live.in, web: www.apexmedicaldevices.com



APEXCARE™
SURGICAL GLOVES

STERILE LATEX
SURGICAL GLOVES



APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES

APEX
MEDICAL DEVICES



STERILE LATEX SURGICAL GLOVES
THE ULTIMATE IN SURGICAL DISPOSABLE PRODUCTS

APEX
MEDICAL DEVICES



CE CERTIFIED COMPANY











CE CERTIFIED COMPANY



100 GLOVES
BY WEIGHT

MADE IN MALAYSIA

Caution - User Should Consider The Circumstances Of Use In Deciding Whether To Remove Any Residual Powder On Gloves After Donning Powder Can Be Removed By Wiping Gloves Thoroughly With A Sterile Wet Sponge, Sterile Wet Towel Or Other Effective Methods."safe Use Of This Gloves By Or On Latex Sensitized Individuals Has Not Established".

Warning - Isolated Cases Pf Allergic Reactions To Rubber Latex Or Powder Have Been Reported. If You Experience A Rection To This Product, Discontinue Use Immediately And Consult Your Physician. This Product Contains Natural Rubber Latex Which May Cause Allergic Reactions In Some Individuals."other Components Used On Making Glove May Cause Allergic Reactions In Some Users".



LOT Batch No. A-2205

Mfg. Date : MAY-2022

Exp. Date : APRIL-2025

Max M.R.P.: ₹ 1499/-

Net 12.61 Gram



100
GLOVES
By Weight

L

POWDER FREE



BLUE



CLEAR

3PTM

Delivering Medical Solutions Worldwide

**VINYL
EXAMINATION
GLOVES**



Scanned with CamScanner

• Non-Sterile • Ambidextrous • Single Use Only



ECONOMY DISPENSER PACK 100 IN BAG

DISPOSABLE PLASTIC GLOVES

- PROTECT HANDS FROM MESSY OILS
- PROTECT FOODS FROM GERMS PERC
- THIN FOR FINGER

NO TASTE
ON OILS

SIZE





APEXCARETM

Disposable Plastic Examination
GLOVES



CAUTION :

Do not use if package is
opened or damaged.

QTY. : 1 PCS.



LOT BATCH NO. : DPE01




MFG./STR. DATE : 08/2020



EXP. DATE : 07/2025

MAX. R.P. ₹. : 10 = 00
(Incl. of all Taxes)

 **Mfr.:**

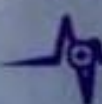

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park,
Opp. State Bank of India, Phase-1, GIDC
Valva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333,
Email : apexmedical@live.in,
web: www.apexmedicaldevices.com

CE

ISO 13485 : 2016

ISO 9001 : 2008

TEAR HERE

APEXCARE™ 1 Unit
NASAL CANNULA

FOR SINGLE PATIENT USE.

Directions for use: Attach oxygen supply tubing to gas source and set the oxygen to the prescribed flow. Check for the gas flow through any device. Insert the nasal tips into the nostrils with two plastic tubes over the ears and under the chin. Then gently adjust the plastic slide until the cannula is secure.

Caution : Oxygen concentration is Affected by the patients tidal volume.

Oxygen liter/min.	% of oxygen
2	28
4	36

SIZE	
ADULT	✓
PEADIATRIC	
INFANT	



FOR SINGLE USE ONLY

LOT : NC-0 5
EXP : 0 5 / 2021
EXP : 0 4 / 2024
MRP : 225=00 Rs.

Manufactured By :

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Opp. State Bank of India, Phase-1, GIDC Vatva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333, Email : apexmedical@live.in, domestic@apexmedicaldevices.com web:www.apexmedicaldevices.com

Be Sure - Be Safe

APEXCARE™ 1 Unit
NASAL CANNULA

FOR SINGLE PATIENT USE.

Directions for use: Attach oxygen supply tubing to gas source and set the oxygen to the prescribed flow. Check for the gas flow through any device. Insert the nasal tips into the nostrils with two plastic tubes over the ears and under the chin. Then gently adjust the plastic slide until the cannula is secure.

Caution : Oxygen concentration is Affected by the patients tidal volume.

Oxygen liter/min.	% of oxygen
2	28
4	36

SIZE	
ADULT	✓
PEADIATRIC	
INFANT	



FOR SINGLE USE ONLY

LOT : NC-0 5
EXP : 0 5 / 2021
EXP : 0 4 / 2024
MRP : 225=00 Rs.

Manufactured By :

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Opp. State Bank of India, Phase-1, GIDC Vatva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333, Email : apexmedical@live.in, domestic@apexmedicaldevices.com web:www.apexmedicaldevices.com

Be Sure - Be Safe

Chlorhexidine/
Ethyl Alcohol
Antiseptic Solution



CHG
HANDRUB



Rub-in
Hand Disinfectant
Kill 99.9 % Germs



FOR CERTIFIED COMPANY
500ml



FOR EXTERNAL USE

Easy to use skin antiseptic
effective against bacteria,
Fungi and viruses.

Precautions / Contraindications

- Do not mix with any other
Chemicals. Keep out of reach of children.
- Keep out of reach of children.
- Contraindicated in persons with
hypersensitivity to Chlorhexidine.

Storage: - Store in cool, dry place
and protect from light.

Topical Antiseptic

Mfg. Lic. No.:
Batch No.:
Mfg. Date:
Exp. Date:
M.R.P. Rs.:
(Incl. of all taxes)



Composition :

Each ml Contains:

Chlorhexidine Gluconate Solution I.P. 2.5% w/v
equivalent to Chlorhexidine Gluconate 0.5% w/v
Ethyl Alcohol IP (with Denaturant) 70% v/v
Colour : Carmoisine

Skin Protecting Substances / Emollients

Application / Directions for use :
Use undiluted

*Use 3ml of (2pushes of dispenser)
on clean, dry hands for 30 second for
Hygienic Hand Disinfection.
Rub 9ml (6 pushes of dispenser) up to
elbows for surgical Hand Disinfection
for 1.5 to 3 minutes

For Rapid Hand Antisepsis :

Dispense sufficient Chlorcare Handrub to
wet hands and finger nails. Rub vigorously
over all surfaces until dry.

For Skin Antisepsis :

Dispense sufficient Chlorcare. Hand rub to
wet area to be treated. Rub vigorously over all
surface until dry.

Chlorhexidine/ Ethyl Alcohol Antiseptic Solution



**CHG
HANDRUB**



CE CERTIFIED COMPANY

#SDB###



**Rub-in
Hand Disinfectant
Kill 99.9 % Germs**



FOR EXTERNAL USE ONLY

Easy to use skin compatible hand rub,
effective against bacteria including bacillus,
Fungi and viruses

Precautions / Contraindications :

- Do not mix with detergents of other
Chemicals. Keep out of eyes and ears.
- Keep out of reach of children.
- Contraindicated in case of known
hypersensitivity to Chlorhexidine.

Storage : Shall be stored below 25° C
and protect from light.

Topical Antiseptic Solution

Mfg. Lic. No.: G/25/2124

Batch No.:

Mfg. Date: AXC-002

Exp. Date: 07-2020

M.R.P. Rs.: 06-2022

(Incl. of all taxes)

5 LTR.

Manufactured for:
APEX MEDICAL DEVICES
Plot - 108, Shree Ujjain Industrial Park,
Opp. S.B.I. Bank, Phase 1, G.I.D.C. Vadod,
Ahmedabad-382445 Gujarat, India.
E-mail : apexmedicaldevices@gmail.com
domestic@apexmedicaldevices.com
Web : www.apexmedicaldevices.com
Mfg. by : Apex Pharmaceuticals Pvt. Ltd. (PVT.)
Haty Chemicals Approach Road, Vadod



Apexcare
Topical Antiseptic Solution
100 ML
99.9% GERM FREE
Rub-in Hand Disinfectant
Liquid hand disinfectant, effective against Bacteria, Fungi & Viruses
ISO CERTIFIED COMPANY



APEX CARE™

Topical
Antiseptic
Solution

HANDRUB

**99.9%
GERM
FREE**

500 ML



**Rub-in
Hand Disinfectant**
Liquid hand disinfectant,
effective against
Bacteria, Fungi & Viruses


ISO CERTIFIED COMPANY



PRECAUTIONS

- Keep away from
- Avoid direct
- Store in cool

DIRECTION FOR USE

Hygienic Hand Disinfection
Approx. 2-3 ml of the solution
rubbed for 20 seconds over
well over clean hands and
nail grooves.

Do not use with

Mfg. Lic. No.:
Batch No.:
Mfg. Date:
Exp. Date:
M.R.P. Rs.:
(Incl. of all taxes)



APEXCARE™

100
BLADES

CARBON STEEL STERILIZED SURGICAL BLADES

Individual Foil Peel Packs
V.C.I Protected
Sterilized by
Gamma Radiation 2.5M-Rads



CE CERTIFIED COMPANY

Skalpell-Klingen
Lames de Bistouri
Hojas de Bisturi

MADE IN INDIA

Fitment No:3



10



11



11P



12



15

Fitment No:4



20



21



22



23



24



APEX CARE
COTTON CREPE BANDAGE
Width : 6 CM x 4 mtrs.
STRETCHED LENGTH
Flesh Coloured Bandage
Washing Renews Elasticity
WITH FAST EDGES
Wig. for :
APEX MEDICAL DEVICES
Plot No. 105, Shyam Ujjawal Ind. Park, Opp. State Bank of India,
G.D.C. Vatva, Ahmedabad-382445, Gujarat, India
Customer Care : +91 79-29701333, Email : apexmedical@gmail.com
StarCott Health Care, Dist. Gandhinagar, Gujarat

APEX CARE
COTTON CREPE BANDAGE
Width : 8 CM x 4 mtrs.
STRETCHED LENGTH
Flesh Coloured Bandage
Washing Renews Elasticity
WITH FAST EDGES
Wig. for :
APEX MEDICAL DEVICES
Plot No. 105, Shyam Ujjawal Ind. Park, Opp. State Bank of India,
G.D.C. Vatva, Ahmedabad-382445, Gujarat, India
Customer Care : +91 79-29701333, Email : apexmedical@gmail.com
StarCott Health Care, Dist. Gandhinagar, Gujarat

APEX CARE
COTTON CREPE BANDAGE
Width : 10 CM x 4 mtrs.
STRETCHED LENGTH
Flesh Coloured Bandage
Washing Renews Elasticity
WITH FAST EDGES
Wig. for :
APEX MEDICAL DEVICES
Plot No. 105, Shyam Ujjawal Ind. Park, Opp. State Bank of India,
G.D.C. Vatva, Ahmedabad-382445, Gujarat, India
Customer Care : +91 79-29701333, Email : apexmedical@gmail.com
StarCott Health Care, Dist. Gandhinagar, Gujarat

APEX CARE
COTTON CREPE BANDAGE
Width : 15 CM x 4 mtrs.
STRETCHED LENGTH
Flesh Coloured Bandage
Washing Renews Elasticity
WITH FAST EDGES
Wig. for :
APEX MEDICAL DEVICES
Plot No. 105, Shyam Ujjawal Ind. Park, Opp. State Bank of India,
G.D.C. Vatva, Ahmedabad-382445, Gujarat, India
Customer Care : +91 79-29701333, Email : apexmedical@gmail.com
StarCott Health Care, Dist. Gandhinagar, Gujarat





APEXCARETM
ELASTIC ADHESIVE BANDAGE B.P.

Width : 10 cm x 4 / 6 mtrs.
STRETCHED LENGTH

Mfg. for :
APEX MEDICAL DEVICES

Plot No.10/B, Shyam Ujjawal Ind. Park, Opp. State Bank of India,
Phase-1, GIDC-Vatva, Ahmedabad-382445, Gujarat, India
Customer Care : +91 79-29701333, Email : apexmedical@live.in
Mfg. By : SteriCott Health Care, Dist: Gandhinagar. (Gujarat)
Lic. No. : G/25/2104



ADULT ☐

Breathe.....S Made Easy

PAEDIATRIC ☐

OxyCare

OxyCare

Multi Flow Venturi Oxygen Mask

- Single use only
- Store in cool and dark place
- Under guidance of qualified medical doctor
- Device to be used for short-term only
- This device to be sold by or on order of Registered Medical practitioners only.
- Discard of device liable to infection may cause to serious illness or death.

LOT : 2006003

Mfg. Dt. : JUNE 2020

Exp. Dt. : MAY 2024

650/- Rs. Per Piece



DEHP



APEX MEDICAL DEVICES

10/B, SHYAM UJJAWAL IND. PARK,
PHASE-1, G.I.D.C. VATVA, AHMEDABAD-382445 (INDIA)

Customer Care No. +91-297011333

Email : apexmedical@live.in







3 BALL SPIROMETER



CE CERTIFIED COMPANY






1 PCS.






FITTING THE MASK

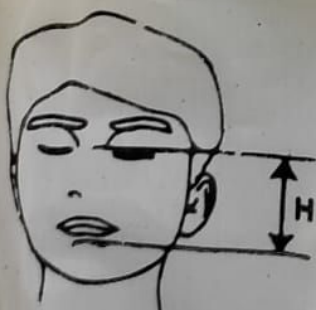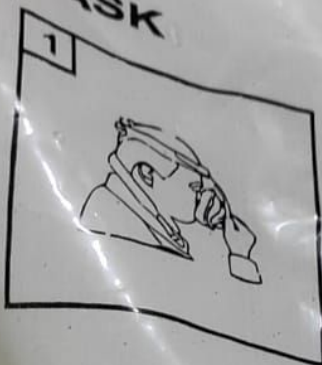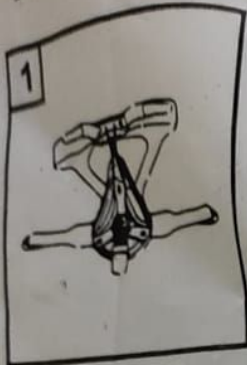

H(mm)	Size
$80 \leq H \leq 90$	S
$90 < H \leq 100$	M
$100 < H \leq 115$	L

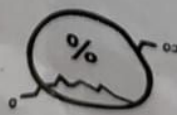



C-PAP/ BIPAP FULL FACE MASK NON VENTED
BATCH No. : 21F035
MFG DATE : 06/2021
EXP DATE : 05/2024
SIZE : SMALL
MRP : 6000/-
(Incl. of All Taxes)
Manufactured By :
APEX MEDICAL DEVICES
Plot No 10/8, Sanyam Ujwal Ind Park,
G.I.D.C, Vaha, Ahmedabad - 382445, Gujarat, India
Customer Care No +91-79-29701333, +91-98246 99959
Email: apexmedical@live.in, domestic@apexmedicaldevices.com
Web: www.apexmedicaldevices.com



FITTING THE MASK

APEX CARE YANKAUR SUCTION SET

INSTRUCTIONS TO USE

- Connect yankaur suction set to suction pump.
- Close vacuum control window by placing thumb on it.
- Insert tip of handle in operative body part and let it remove secretion by sucking.
- Vacuum can be controlled on regular intervals by vacuum control window.

CE

ISO 13485 : 2018
ISO 9001: 2008

CAUTION


- Do not Reuse.
- Do not use if package is open or damaged.
- Do not store at extreme temperature and humidity.
- Discard after use.
- Use under guidance of qualified medical practitioner.

E.T.O STERILISED | NON TOXIC | PYROGEN FREE | DISPOSABLE

Mfr.:

APEX MEDICAL DEVICES

Plot No. 10/B, Shyam Ujjawal Ind Park,
Opp. Late Bank of India, Phase-1, GIDC-Vatva,
Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333.

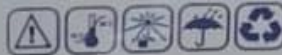

2

STERILE EO

Qty : 1 Unit.

LOT / BATCH NO.
MFG. / EST. DATE
EXP. DATE
MAY 19 2025
(Valid for all States)

YEAR HERE

APEXCARE™

NEBULIZER FACE MASK

- DISPOSABLE
- FOR SINGLE PATIENT USE
- NON INTENDED FOR REPROCESSING

**ACCESSORY FOR COMPRESSED
AIR BREATHING APPARATUS**

CE

ISO 13485 : 2016
ISO 9001: 2008

- Discard after single use.
- Jeter après utilisation.
- Desechard duspues de utizar.
- Zum einmalgebrauch.

ADULT ☐

CHILD ☐

STORE IN DRY AND COOL PLACE

Manufactured By:

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Opp. State Bank of India,
Phase-1, GIDC Vathe, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333, Email: apexmedical@live.in,
domestic@apexmedicaldevices.com web: www.apexmedicaldevices.com

Lot No. :
Mfg. Date :
Exp. Date :
M.R.P
(Incl. of all Taxes)



FOR SINGLE USE ONLY
Quantity : 1 Unit

Be Sure - Be Safe

TEAR HERE

APEXCARETM

OXYGEN FACE MASK

- DISPOSABLE
- FOR SINGLE PATIENT USE
- NON INTENDED FOR REPROCESSING

**ACCESSORY FOR COMPRESSED
AIR BREATHING APPARATUS**

- Discard after single use.
- Jeter après utilisation.
- Desechard duspues de uñizar.
- Zum einmalgebrauch.

STORE IN DRY AND COOL PLACE

CE

ISO 13485 : 2016
ISO 9001: 2008

ADULT ☐

CHILD ☐

Manufactured By :

APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Opp. State Bank of India,
Phase-1, GIDC Vatva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333, Email : apexmedical@live.in,
domestic@apexmedicaldevices.com web:www.apexmedicaldevices.com

Lot No. :

Mfg.Date :

Exp.Date :

M.R.P ₹. :
(Incl. of all Taxes)



FOR SINGLE USE ONLY
Quantity : 1 Unit

Be Sure - Be Safe

Tear Here

APEXCARETM

NASAL CANNULA

1 Unit

FOR SINGLE PATIENT USE.

Directions for use: Attach oxygen supply tubing to gas source and set the oxygen to the prescribed flow. Check for the gas flow through any device. Insert the nasal tips into the nostrils with two plastic tubes over the ears and under the chin. Then gently adjust the plastic slide until the cannula is secure.

Caution : Oxygen concentration is Affected by the patients tidal volume.

Oxygen liter/min.	% of oxygen
2	28
4	36

SIZE	
ADULT	✓
PEADIATRIC	
INFANT	



FOR SINGLE USE ONLY

LOT : NC-0 5
 : 0 5 / 2021
 : 0 4 / 2024
MRP. : 225=00 Rs.
(Incl. All taxes)

Manufactured By :



APEX MEDICAL DEVICES

Plot No 10/B, Shyam Ujjawal Ind Park, Opp. State Bank of India,
Phase-1, GIDC Vatva, Ahmedabad-382445, Gujarat, India.
Customer Care: +91-79-29701333, Email : apexmedical@live.in,
domestic@apexmedicaldevices.com web:www.apexmedicaldevices.com

Be Sure - Be Safe





SURGICAL DISPOSABLE FACE MASK



CE CERTIFIED COMPANY



50 Pcs

Disposable Face Mask

Surgical Face Mask - 3 Layers Structure



Ear Strap on both ears



Adjust the bridge of nose



Unfold the folding part



SURGICAL DISPOSABLE FACE MASK

Matters Needing Attention

- This Mask should be replaced in time, not recommended for long-term use.
- It is recommended to stop using the product in case of any adverse reaction during wearing.
- This product cannot be washed. Please make sure it is used within the validity period.
- Store in a dry and ventilated place away from fire sources and inflammables.



CE CERTIFIED COMPANY



50 Pcs





Pack Of 100







FILTRATION RESPIRATOR

DUSTPROOF / ANTI-FOG

PROTECTION

Disposable Foldable mask

N-95 Standard BFE $\geq$ 95%



**5
Pcs.**



N95 MASK

FILTRATION RESPIRATOR

DUSTPROOF / ANTI-FOG

PROTECTION

Disposable Foldable mask

N-95 Standard BFE $\geq$ 95%



**5
Pcs.**







Face Shield

High-visibility
Ultra-clear
Easy to use
All-weather
Easy to clean
All-weather
Easy to use





HIV-KIT



CE CERTIFIED COMPANY

MFR:
APEX MEDICAL DEVICES

Phase - 1, GIDC - Vatva,
Ahmedabad - 382445, Gujarat, India.
Customer Care : +91-79-29701333,
Website : www.apexmedicaldevices.com
Email : apexmedical@live.in



Hand Gloves



3 PLY Mask



Face Shield



Disposable Bag



Shoe Cover

Personal Protection Equipment Kit



APEXCARE™
PPE KIT

CONTENTS:

- | | |
|------------------------------------|-----------|
| 1) Coverall Suits | : 1 Piece |
| 2) 3 Layer Face Mask / N 95 | : 1 Piece |
| 3) Sterile Gloves / Nitrile Gloves | : 1 Pair |
| 4) Safety Goggles / Face Shield | : 1 Piece |
| 5) Disposable Waste Bag | : 1 Piece |
| 6) Shoe Cover | : 1 Pair |

Lot No. : PPE 01
MFG. DT. : May-2020
Exp. DT. : April-2023
Max Retail Price Rs. : 2599.00
(Incl. All taxes)

PACKED & MARKED BY:
APEX MEDICAL DEVICES
Phase - 1, GIDC - Vatva,
Ahmedabad - 382445, Gujarat, India.
Customer Care : +91-79-28701333,
Website : www.apexmedicaldevices.com
Email : apexmedical@live.in

